



PurMop BLACK® IF9-S

Cleanroom Mop for ICT (Isolator
Cleaning Tool) | Performance Line

Sterile | Fillpack (ready to fill with fluid)

Polyester-Microfiber

Item no. 2112050

STERILE

ISO	3	4	5	6	7	8
GMP			A/B		C	D

The IF9-S is a sterile and chemical resistant pouch containing 2 cleanroom mops for ICTs* – ready for saturation. The pouch can easily be filled with the fluid of choice (common bactericides, sporicides and detergents used in cleanrooms). Thanks to their high-purity knitted microfiber fabric, the mops ensure excellent dirt pick-up, an even liquid release and very low particle emissions. The mops have been developed to match the PurMop BLACK® carbon fiber isolator cleaning tools. Special material flaps ("wings") allow easy corner cleaning.

- Chemical-resistant pouch with screw cap: ready for saturation with the fluid of choice
- "BLACK21" – new mop size suitable for PurMop BLACK® isolator cleaning tools
- Microfiber with high-low structure for optimized cleaning performance
- With "wings" for simplified corner cleaning

Application workbenches / work surfaces

Industry sectors pharma / medical /
pharmacy / laboratory / mechanics

*mop model:
PurMop BLACK® ID9-S (item no. 2112034)
Request data sheet for detailed information.

APPLICATION, USAGE

Recommended for cleanroom cleaning and disinfection in isolators, workbenches and critical production environments as well as difficult to reach areas on machines and equipment.

SIZE, WEIGHT, COLOR

8.27 x 4.33 " (±0.2") | 0.8 oz (± 10%)
white with orange loop

FRONT MATERIAL

100% polyester-microfiber

BACKING MATERIAL

100% polyester

EDGES

100% polyester tape (bound)

FILLING QUANTITY

3.38 US fl.oz.

SURFACE COVERAGE

Up to 129 ft² (per fillpack)

CHEMICAL COMPATIBILITY

Against acids: resistant
Against alkalis: resistant

ABRASION RESISTANCE

≥ 20.000 movements (9 kPa)

STERILITY

Gamma irradiated, SAL 10⁻⁶,
validated according to ISO 11137

SHELF LIFE

2 years

RE-USE, RECYCLING

Mops recyclable

STORAGE

General storage conditions apply, protect from harmful influences (fluids, sun rays or UV light, contamination, pressure), recommended storage temperature 41–104 °F

PACKAGING MATERIAL

Plastic pouch, white: PET/PE
Screw cap, white: PE

PACKING

2 mops / pack
30 packs / case [60 mops]

CASE SIZE (L/W/H)

CASE WEIGHT
23.23 x 15.75 x 6.30" | 7.3 lbs

PALLET

48 cases / pallet
(47.24 x 33.5 x 83.86")

COUNTRY OF ORIGIN

Socialist Republic of Vietnam

COMPATIBLE ACCESSORIES

- PurMop BLACK® isolator cleaning tools
- Funnel (diameter 4.72", polypropylene, item no. 2112060)
- Measuring cup (volume 8.45 US fl.oz., polypropylene, printed blue scale, item no. 2112066)

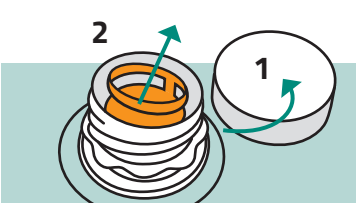
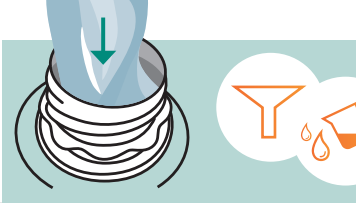
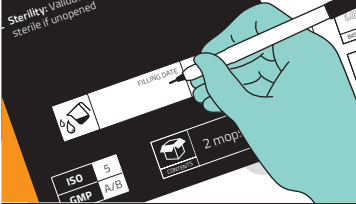
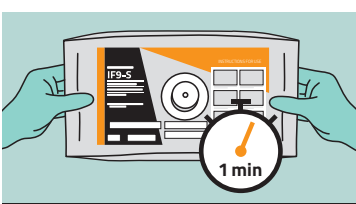
>>



PurMop BLACK®

IF9-s

INSTRUCTIONS FOR USE

NOTES		
The filled pouches must not be stored for longer than 24 hours. Recommended: The compatibility of the chemical formulation of the fluid and the surface to be cleaned or disinfected should be checked before adding the liquid to the dry mops.		
STEP 1	STEP 2	
		Unscrew the lid (1) and remove sealing (2). Add 3.38 US fl.oz. of fluid and close the lid.
STEP 3	STEP 4	
		Note the fluid used and the filling date on the label. Wait for 1 minute
STEP 5	STEP 6	
		Open pouch by tearing from edge Pick up mop with PurMop BLACK® isolator cleaning tool